WAVEJ®

ETB-129237.1

ETB-129237.2

ETB-129237.3

- Elektrische Zahnbürste für Kinder (D)***
Elektrische tandenborstel voor kinderen (NL)
Cepillo de dientes eléctrico para niños (E)
Escova de dentes elétrica para criança (P)
Spazzolino elettrico per bambini (I)
Detská elektrická zubná kefka (SK)
Otroška električna zobna ščetka (SLO)
Elektryczna szczoteczka do zębów dla dzieci (PL)
Periuță de dinți electrică pentru copii (RO)
Dječja električna četkica za zube (HR)
Elektrický zubní kartáček pro děti (CZ)
Електрическа четка за зъби за деца (BG)
Gyermek elektromos fogkefe (H)
Kids electric toothbrush (EN)

Inhalt – Inhoud – Contenido – Índice – Contenuto – Obsah – Vsebina – Treść –
Cuprins – Sadržaj – Obsah – Съдържание – Tartalom – Content

Bedienungsanleitung – German (D)	- 2 -
Gebruiksaanwijzing – Netherlands (NL)	- 8 -
Manual de Instrucciones – Spanish (E).....	- 14 -
Manual de instruções – Portuguese (P)	- 20 -
Manuale d'istruzioni – Italian (I).....	- 26 -
Návod na obsluhu – Slovak (SK)	- 32 -
Navodila za uporabo – Slovenian (SLO).....	- 38 -
Instrukcja obsługi – Polish (PL)	- 44 -
Manual de utilizare – Romanian (RO)	- 50 -
Uputa za korištenje – Croatian (HR)	- 56 -
Návod k použití – Czech (CZ)	- 62 -
Ръководство с инструкции – Bulgaria (BG)	- 68 -
Kezelési útmutató – Hungary (H)	- 74 -
Instruction manual – English (EN)	- 80 -

Bedienungsanleitung – German (D)

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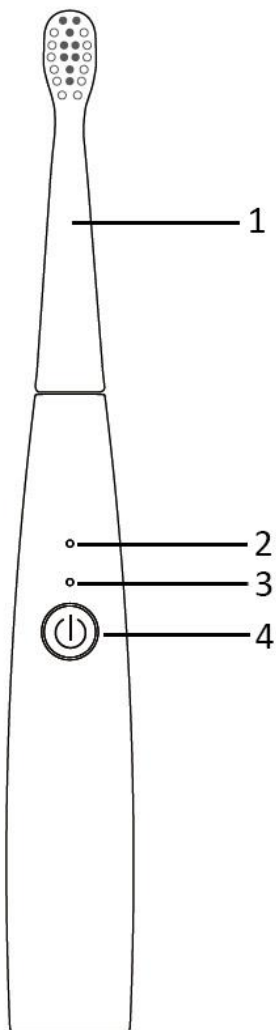
Instruction manual – English (EN)

SAFETY INSTRUCTIONS

Before use make sure to read all of the below instructions in order to avoid injury or damage, and to get the best results from the appliance. Make sure to keep this manual in a safe place. If you give or transfer this appliance to someone else make sure to also include this manual.

In case of damage caused by user failing to follow the instructions in this manual the warranty will be void. The manufacturer/importer accepts no liability for damages caused by failure to follow the manual, a negligent use or use not in accordance with the requirements of this manual.

1. Toothbrushes can be used by children and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Cleaning and user maintenance shall not be made by children without supervision.
2. Children shall not play with the appliance.
3. Check and make sure to use an adapter with a rated output voltage of 5V, a rated output current of 1A, and a waterproof level of IPX4 to ensure the electrical safety. The adapter cannot be used for charging this product if it does not comply with the above requirements.
4. If the appliance drops into water while it is being charged, please firstly disconnect the power source and then remove the appliance from the water.
5. Please charge the appliance in a dry environment far away from water. Do not charge the appliance in a humid environment such as bathroom and others.
6. It is forbidden to use and rinse the toothbrush while it is being charged.
7. Keep your hands dry when handling the power cable for charging.
8. Do not place the appliance and its battery in a microwave oven or on an induction cooker.
9. Do not disassemble or modify the toothbrush before disposing of it for recycling.
10. Do not immerse the toothbrush into hot water.
11. Do not press against tooth or gum with the toothbrush using excessive force. Do not brush the same position for long periods time to avoid any damage.
12. **WARNING:** Make sure that there is no appearance of damage before each use. Do not use the toothbrush if damaged.
13. **WARNING:** Please consult your doctor whether you can use the toothbrush If you have had dentistry surgical operation in past two months or are under being subject to dental treatment.
14. **WARNING:** Oral hemorrhage may occur when using this product at the first time. If the oral hemorrhage lasts for more than a week, please consult your doctor.
15. **WARNING:** Please consult your doctor before using this product if you are using a heart pacemaker or other electronic devices on or in your body.
16. **WARNING:** Please consult your doctor before using this product if you have medical problems.
17. **HAZARD:** Do not heat the toothbrush, throw it into the fire, store, charge or use at a high-temperature environment to avoid overheating, fire or explosion as the handle of this product is equipped with a built-in rechargeable battery.
18. The toothbrush is for personal use, do not share it with other persons.
19. This product is designed for cleaning teeth, gums, and tongue only.
20. The temperature for charging, using, and storing the product shall be between 5°C and 40°C.
21. Please recharge this product in time when the battery is low.
22. Do not use the toothbrush when you are taking a bath or shower.
23. This appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
24. This appliance contains batteries that are non-replaceable.
25. Regarding the information pertaining to the use of the appliance and user maintenance, thanks to refer to the below paragraph of the manual.


PARTS DESCRIPTION

1. Toothbrush head
2. Morning mode indicator light
3. Night mode indicator light
4. Switch button

BEFORE USE**1. Charge your toothbrush with an adapter**

- 1) Turn off the toothbrush before charging.
- 2) Place the toothbrush on a working surface horizontally. Insert one end of the provided USB cable into the DC charging port (a) at the bottom of the toothbrush first and then the other end into an adapter. Connect the adapter with a proper power source.

Note: Check and make sure to use an adapter with a rated output voltage of 5V, a rated output current of 1A, and a waterproof level of IPX4 to ensure the electrical safety. The adapter cannot be used for charging this product if it does not comply with the above requirements.

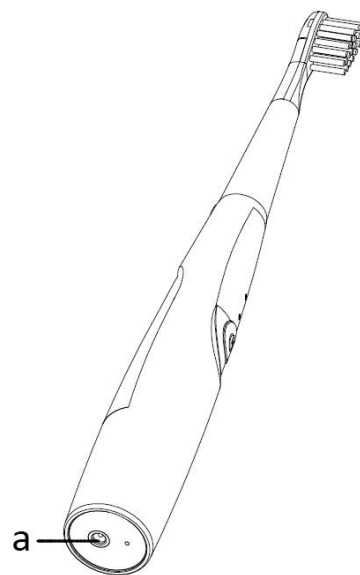- 3) During charging, the morning mode indicator light flashes in white and the night mode indicator light is off. Once the toothbrush is fully charged, the night mode indicator light will illuminate in blue solidly and the morning mode indicator light will be off.

Note: It may take approx. 10 hours to fully charge the toothbrush.

Note: Do not turn on the toothbrush during charging.

Notes:

- When the toothbrush is charged for the first time or if the toothbrush has not been used for a long time, the indicator light may not illuminate for the first 5 minutes during charging. This is caused by the low battery of the toothbrush, which is normal.
- The toothbrush must be charged completely under the following conditions:
 - The toothbrush is used for the first time.
 - The toothbrush has not been used for more than 3 months.
 - The toothbrush vibration becomes weak or low-battery reminding is shown.
- The toothbrush may become hot when it is being charged. This is normal. This will not affect the performance of the toothbrush.


Low-battery reminding function:

During working, the indicator light of the selected mode will flash and then extinguish when the toothbrush is of low battery.


2. Installing / detaching the toothbrush head

Installing: Align the mark A of the toothbrush head with the mark B of the handle. Push the toothbrush head downward to the handle and then rotate clockwise until the mark A is aligned with mark T.

Detaching: Turn the toothbrush head anti-clockwise and then pull it upwards.

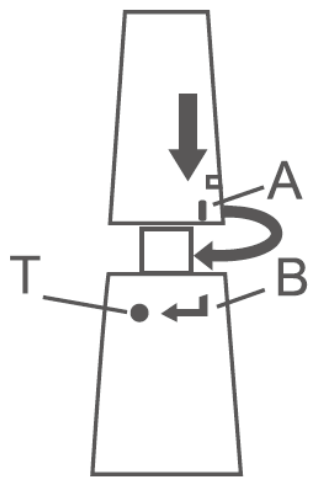
Install

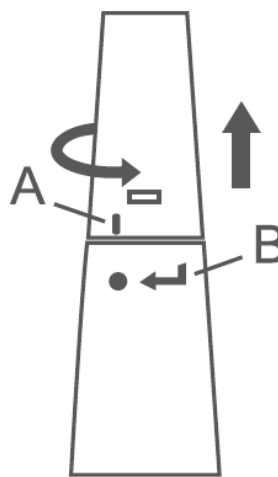Detach

USE

1. Put some toothpaste on the bristles before activating the toothbrush by pressing the switch button.
2. Put the toothbrush into your mouth and then place the toothbrush head between the tooth and the gum.
3. Press the switch button to turn on the toothbrush. The toothbrush will operate as followed.
 - 1) The toothbrush starts running at the morning mode; its indicator light illuminates in white. The toothbrush works at a full power. The vibration frequency and the cleaning efficiency are high. It is recommended to use the morning mode in the morning to clean the dental plaque produced at night.
 - 2) Press the switch button again **WITHIN** 5 seconds to enter the night mode, its indicator light will illuminate in blue. The toothbrush works at a lower power. The vibration frequency is lower. It is recommended to use the night mode before sleep to care for the gum.
Alternatively, press the switch button again **AFTER** 5 seconds, the toothbrush will turn off directly.
 - 3) With no further operations in 5 seconds, the toothbrush will keep running at the current mode and shut off automatically in 2 minutes. The current mode will be saved and the toothbrush will work according to the saved mode when turned on again.
Note: When turned on again and the toothbrush starts running at the night mode, press the switch button again:
 - **WITHIN** 5 seconds, the toothbrush will turn off directly. When turned on again, the toothbrush will work at the morning mode.
 - **AFTER** 5 seconds, the toothbrush will turn off directly. When turned on again, the toothbrush will work at the night mode.
4. When brushing your teeth, gently move the toothbrush head back and forth between the teeth slowly. The toothbrush will pause for a short while approx. every 30 seconds to remind you to brush other sections of your mouth.
Note: There are two rows of teeth in the mouth and each row of teeth have left and right sections, so that there are four sections in the mouth needed to be brushed.
5. When finished cleaning teeth, turn off the toothbrush, rinse the toothbrush head and bristles with running water and then dry them. Wipe the handle with a towel.

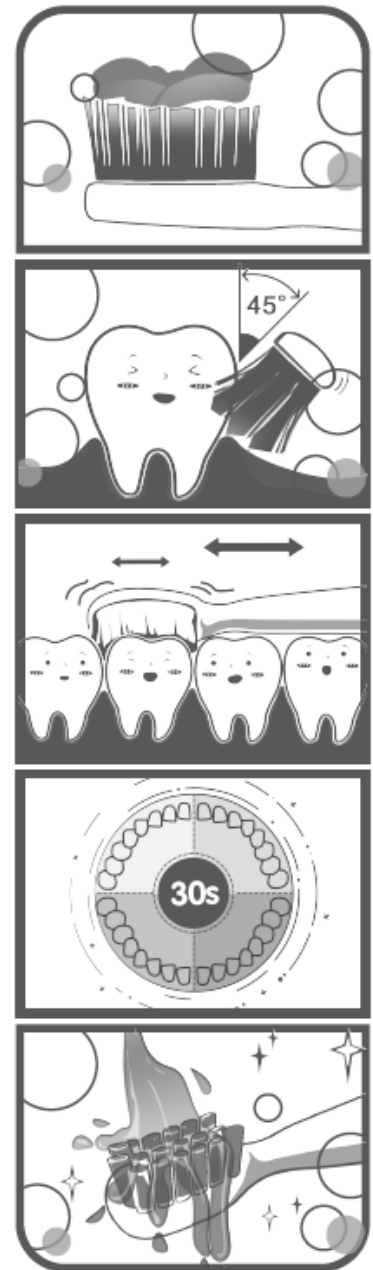CLEANING AND MAINTENANCE

CAUTION: Before cleaning, disconnect the toothbrush from the power source.

Rinse the toothbrush head with clean water at a temperature lower than 45°C and then dry the surface of the toothbrush handle with a towel.

CAUTION: Do not immerse the toothbrush into water. Never use poisonous cleaning agents.

CAUTION: Do not let water get into the charging port on the bottom of the toothbrush when cleaning.

CAUTION: Do not put any part of the product into the dishwasher for cleaning.

CAUTION: Do not use chemical, steel, wooden or abrasive cleaners to clean the toothbrush.

CAUTION: Please do not use a deformed toothbrush head or a toothbrush head with bent bristles.

CAUTION: Replace the toothbrush head every 3 months or sooner if signs of wear appear.

Clean the toothbrush head

1. Rinse the toothbrush head and bristles with running water after each use.
2. Remove the toothbrush head and clean its inside surface with warm water at least once a week.

Clean the handle

1. Remove the toothbrush head and rinse the vibration piece with warm water. Make sure all the residual toothpaste is removed.
2. Wipe the handle surface with a wet cloth and then dry thoroughly.


Storage

If the toothbrush will not be used for a long time, disconnect the power source and clean the toothbrush accordingly. Store the toothbrush in a cool and dry place away from the direct sunlight. Regularly charge the toothbrush every three months to prevent the battery from discharging completely. A fully-discharged battery may not be able to be recharged.

TROUBLESHOOTING

Problem	Cause	Solution
Vibration becomes weak.	1. The toothbrush works with too much pressure against the teeth. 2. Low battery	1. Place the toothbrush gently on the teeth. 2. Please recharge the battery.
The toothbrush only works for a short time after being charged.	The charging time is not enough.	Charge the toothbrush completely for approx. 10 hours.
The toothbrush stops working.	The battery has run out.	Please recharge the battery.

TECHNICAL DATA

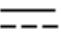Rated input: 5V  0.2A

GUARANTEE AND CUSTOMER SERVICE

Before delivery our devices are subjected to rigorous quality control. If, despite all care, damage has occurred during production or transportation, please return the device to your dealer.

For the purchased device we provide 2 years guarantee, commencing from the day of sale. If you have a defective product, you can directly go back to the point of purchase.

Defects which arise due to improper handling of the device and malfunctions due to interventions and repairs by third parties or the fitting of non-original parts are not covered by this guarantee. Always keep your receipt, without the receipt you can't claim any form of warranty. Damage caused by not following the instruction manual, will lead to a void of warranty, if this results in consequential damages then we will not be liable. Neither can we hold responsible for material damage or personal injury caused by improper use if the instruction manual is not properly executed. Damage to accessories does not mean free replacement of the whole appliance. In such case please contact our service department. Broken glass or breakage of plastic parts is always subject to a charge. Defects to consumables or parts subjected to wearing, as well as cleaning, maintenance or the replacement of said parts are not covered by the warranty and are to be paid.

Electrical and electronic devices/ information for private households

The Directive 2012/19/EU on waste electrical and electronic equipment contains a large number of requirements for handling electrical and electronic equipment. The most important are compiled here.

1. Separate collection of waste electrical and electronic equipment

Electrical and electronic equipment that has become waste is referred to as waste electrical and electronic equipment (WEEE). Owners of WEEE must collect them separately from unsorted municipal waste. WEEE do not belong in household waste, but in special collection and return systems.

2. Batteries, accumulators and lamps

Owners of WEEE must hand over them in complete condition. However, owners of WEEE have to separate old batteries and accumulators which are not enclosed by the WEEE, as well as lamps which can be removed, from the WEEE without destroying it before handing it over to a collection point.

3. Ways to return the WEEE

Owners of WEEE from private households can hand them in free of charge at the designated collection points of the public waste disposal authorities or at the take-back points set up by manufacturers or distributors in terms of the Directive 2012/19/EU.

4. Privacy Protection Notice

WEEE often contain sensitive personal data. This applies especially to information and telecommunications technology devices such as computers and smartphones. In your own interest, please note that each end user is responsible for deleting the data on the WEEE before disposing it.

5. Potential environmental impact

WEEE contain substances that may have a negative impact on the environment and human health if their collection, hand-over, re-use or material recovery is not done in compliance with respective legislation.

6. Your role in WEEE waste treatment

By fulfilling these obligations and especially by fulfilling the obligation to collect WEEE separately, not to dispose it together with unsorted municipal waste and to hand it over to at collection points, you as the end user contribute to the re-use and material recovery of WEEE.

7. Meaning of the crossed-out dustbin symbol

The crossed-out dustbin, which is regularly displayed on electrical and electronic equipment, indicates that the respective device must be collected separately from unsorted municipal waste at the end of its life cycle.

BATTERY MUST BE RECYCLED OR DISPOSED OF PROPERLY. DO NOT OPEN. DO NOT DISPOSE OR IN FIRE OR SHORT CIRCUIT.

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Emerio B.V.