emerio®

SMG-124653

Massage Gun (EN)
Massagepistol (SE)

CE

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Instruction manual – English	- 2 -
Bruksanvisning – Swedish.....	- 7 -

SAFETY INSTRUCTIONS

Before use make sure to read all of the below instructions in order to avoid injury or damage, and to get the best results from the appliance. Make sure to keep this manual in a safe place. If you give or transfer this appliance to someone else make sure to also include this manual.

In case of damage caused by user failing to follow the instructions in this manual the warranty will be void. The manufacturer/importer accepts no liability for damages caused by failure to follow the manual, a negligent use or use not in accordance with the requirements of this manual.


1.  Read operator's manual.
2. This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.
3.  Detachable supply unit <XH2520-0450wG>
WARNING: For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.
4. This appliance contains batteries that are non-replaceable.
5. The appliance is only to be used with the power supply unit provided with the appliance.

Forbiddance:

1. It is forbidden for minors to use. If the muscles have been damaged, bruises, and there are signs of pain or discomfort, the usage of this product is prohibited. Please consult your doctor before using this device.

2. It is forbidden to use on non-dry, non-clean skin surface. You can use it across clothes, gently press and move on the skin, each part for about 60 seconds.
3. It is forbidden to use on areas above the neck and with less muscle coverage, on bones, and only use this product on the soft tissues of the body without pain or discomfort.
4. It is forbidden to disassemble this product in any way.
5. It is forbidden to place any objects in the product ventilation holes.
6. It is forbidden to immerse the product in water, throw the product into the fire or close to the fire source. It is forbidden to drop the product intentionally.
7. It is forbidden to use other chargers for charging. You must use the charger provided with the appliance.
8. It is forbidden to keep fingers, hair or other parts of the body near the shaft and back of the massage head, because it may be caught or burned. While installing the massage head, do not touch the installation part.

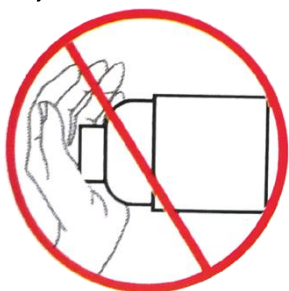Do not touch the device when massage gun is working.

Warning:

1. Keep away from liquids or heat sources and keep the vents free of dust.
2. When charging is complete, unplug the charger.
3. To use the product during charging is prohibited.
4. Please use this product according to the instructions.
5. Before each use, carefully check whether the product can be used normally according to the operation method of this manual.
6. Avoid exposing the machine to high temperature.

Caution:

This product cannot be used in any of the following situations without the prior consent of the doctor.

Pregnancy, diabetic complications (such as neuropathy or retinal damage), wearing a pacemaker, recent surgery, epilepsy or migraine, disc herniation, spondylolisthesis, vertebral detachment, or spinal joint rigidity, joint replacement surgery or IUD, metal needles or plates, and any problems related to your physical health. These contraindications do not mean that you cannot use this product. We recommend you consult with a doctor before use. The ongoing research is about the impact of vibration training on specific medical diseases. This is likely to reduce the above contraindications. Practical experience shows that in many cases, combining vibration training with a treatment plan is also a treatment method. This must be done on the advice of a doctor, expert or physiotherapist, or accompanied by them.

PARTS DESCRIPTION

1. Massage head
2. Speed display lights
3. Speed control button
4. Ventilation holes
5. Handle
6. Battery indicator light
7. Charging port
8. Power switch on/off
9. Small round head (Mainly used for large muscles)
10. U-shaped head (Mainly used for neck and spine)
11. Tower-shaped head (Mainly used for soft tissues and sensitive muscle groups)
12. Flat head (Mainly used for various parts of the body)
13. Cylinder head (Mainly used for deep muscle groups)
14. Round head (Mainly used for various parts of the body)
15. Charger

Massage head installation method:

While installing/pulling, the appliance should be in the standby/shutdown state.

- Push the massage head into the appliance as far as possible to install it firmly.
- Slightly rotate the massage head in any direction to pull it out from the appliance.

**USAGE**

1. This massage gun is suitable for physical massage before or after exercise.
2. Turn on the power switch; the green indicator light is on and the appliance is in standby mode.
3. Press the speed control button to enter the first speed level; the appliance begins to run. Each time the button is pressed, the speed level is increased by one level. Press the speed control button at the fifth speed level, the appliance enters the standby mode. Five speed display lights correspond to the five speed levels.
4. Turn off the power switch in any speed level, the appliance will stop immediately.
5. After use, please turn off the power switch. Do not keep it in standby mode for a long time.

Note: This product will switch to the standby mode after running for 10 minutes. To continue using please press the speed control button again.

CHARGING

1. Before initial use, charge for more than 5 hours to fully charge the battery.
2. To start charging, please connect the charger to the battery charging port and the power outlet.
3. The indicator light shows red for charging.
4. When the indicator light turns to green, it means the power is sufficient.
5. The battery can be charged at any time and at any power level. The appliance cannot be turned on during charging.
6. It needs to be charged when the battery is exhausted and the red light is displayed.

MAINTENANCE AND STORAGE

1. Clean the appliance with a slightly damp towel, and then dry it with a soft cloth.
2. When storing or travelling, put the appliance in the original packaging.
3. The appliance should be stored in a dry, well-ventilated place without corrosive gas.
4. After a long period of non-use, it should be cleaned. Remember to recharge the appliance at least every six months.
5. This appliance is suitable for transportation in cars, trains, ships and airplanes.
6. The appliance should be fixed during transportation to prevent violent vibration, collision, rolling and falling, and should be upward and waterproof packed.

TECHNICAL DATA

Voltage: 24V DC

Battery capacity: 2000mAh

Charger input: 100-240V ~ 50/60Hz, 0.5A

Charger output: 25.2V $\overline{\overline{=}}$ 0.45A

GUARANTEE AND CUSTOMER SERVICE

Before delivery our devices are subjected to rigorous quality control. If, despite all care, damage has occurred during production or transportation, please return the device to your dealer. In addition to statutory legal rights, the purchaser has an option to claim under the terms of the following guarantee:

For the purchased device we provide 2 years guarantee, commencing from the day of sale. If you have a defective product, you can directly go back to the point of purchase.

Defects which arise due to improper handling of the device and malfunctions due to interventions and repairs by third parties or the fitting of non-original parts are not covered by this guarantee. Always keep your receipt, without the receipt you can't claim any form of warranty. Damage caused by not following the instruction manual, will lead to a void of warranty, if this results in consequential damages then we will not be liable. Neither can we hold responsible for material damage or personal injury caused by improper use if the instruction manual is not properly executed. Damage to accessories does not mean free replacement of the whole appliance. In such case please contact our service department. Broken glass or breakage of plastic parts is always subject to a charge. Defects to consumables or parts subjected to wearing, as well as cleaning, maintenance or the replacement of said parts are not covered by the warranty and are to be paid.

ENVIRONMENT FRIENDLY DISPOSAL

Recycling – European Directive 2012/19/EU

This marking indicates that this product should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle

it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

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