emerio[®]

MT-105867.2

Mouth thermometer (EN)

Munntermometer (SE)

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CE

WARNING:


1. Choking hazard: thermometer cap and battery may be fatal if swallowed. Do not allow children to use this device without parental supervision.
2. Do not use thermometer in ear. Designed use is for oral, rectal, and armpit (axilla) readings only.
3. Do not place thermometer battery near extreme heat as it may explode.
4. Remove battery from the device when not in operation for a long time.
5. Note: use of the probe cover may result in a 0.1°C (0.2°F) discrepancy from actual temperature.
6. It is recommended the performance should be checked every two years.

Consult Accompanying Documents

Please read carefully before using

This digital thermometer provides a quick and highly accurate reading of an individual's body temperature.

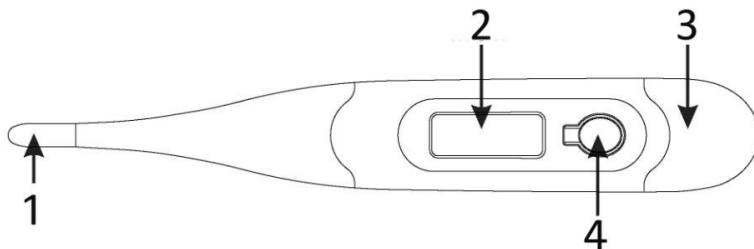PRODUCT ILLUSTRATION

1. Probe
2. LCD
3. Battery Cap
4. On/Off Button

PRECAUTION

Do not bend or drop thermometer. It is not shockproof.

Do not boil probe to disinfect.

Do not store in direct sunlight or in areas of extreme heat.

Do not disassemble. See BATTERY REPLACEMENT to replace battery.

Avoid strong electromagnetic interference such as microwave ovens and cell phones.

SPECIFICATIONS

Range:	32.0°C—42.9°C (90.0°F -109.9°F) chosen by manufacturer (°C/°F)
Accuracy:	±0.1°C, 35.5°C—42.0°C (±0.2 °F, 95.9°F -107.6°F)
	±0.2°C under 35.5°C or over 42.0°C
	(±0.4°F under 95.9°F or over 107.6°F)
	at standard room temperature of 18°C to 28°C (64.4°F to 82.4°F)
Display:	Liquid crystal display, 3 1/2 digits
Memory:	For storing the last measured value
Battery:	One 1.5V DC. button battery (LR41 or SR41, UCC 392)
Battery life:	Approx 200 hours
Dimension:	12.3cm x 1.9cm x 1.1cm (L x W x H)
Weight:	Approx. 11 grams including battery

Using environment:	Temperature: 5°C~40°C (41°F~104°F) Relative humidity: ≤85%RH
Storage and transportation condition:	Temperature: -10°C ~55°C (14°F~131°F) Relative humidity: ≤95%RH
Waterproof level	IP27
Classification:	Type BF 

DIRECTIONS

- Wipe or dip thermometer probe with isopropyl rubbing alcohol to disinfect.
- Press the On/Off button next to LCD display. A tone will sound as the screen shows, followed by last record temperature. After showing the self-test temperature, the thermometer is now in the testing mode.
Note: If the measured temperature is less than 32°C or 90°F, the LCD will display Lo, if the measured temperature is higher than 42.9°C or 109.9°F, the LCD will display Hi.
- Position thermometer in desired location (mouth, rectum, or armpit.)
 - Oral Use: Place thermometer under tongue as indicated by position shown in Figure 2. Close your mouth and breathe evenly through the nose to prevent the measurement from being influenced by inhaled/exhaled air. Normal temperature between 35.7°C and 37.3°C (96.3°F and 99.1°F)
 - Rectal use: Lubricate silver probe tip with petroleum jelly for easy insertion. Gently insert sensor approximately 1cm (less than 1/2") into rectum. Normal temperature between 36.2°C and 37.7°C (97.2°F and 99.9°F).
 - Armpit use: Wipe armpit dry. Place probe in armpit and keep arm pressed firmly. From a medical viewpoint, this method will always provide inaccurate readings, and should not be used if precise measurements are required. Normal temperature is between 35.2°C and 36.7°C (95.4°F and 98.1°F).
- The degree sign flashes throughout the testing process. When flashing stops an alarm will beep for approximately 10 seconds. The measured reading will appear on the LCD simultaneously. The displayed temperature will not change when thermometer is removed from its testing position.
Note: normally the buzzes are "Bi-Bi-Bi-Bi-"; Alarm beeps more rapidly when temperature reaches 37.8°C (100°F) or higher, and the buzzes are "Bi-Bi-Bi---- Bi-Bi-Bi--- Bi-Bi-Bi".
- To prolong battery life, press the On/Off button to turn unit off after testing is completed. If no action is taken, the unit will automatically shut off after around 10 minutes.

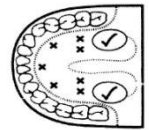
Figure 2

BATTERY REPLACEMENT

- Replace battery when "" appears in the lower right corner of LCD display.
- Pull battery cover off as shown in Figure 3.
- Gently pull out plastic circuit board with battery chamber approximately 1 cm (slightly less than 1/2" (See Figure 4).
- Use pointed object such as a pen to remove old battery. Discard battery lawfully. Replace with new 1.5 DC button type LR41 or SR41, UCC392, or equivalent. Be sure battery is installed with "+" polarity facing up. (See Figure 5)
- Slide battery chamber back into place and attach the cover.

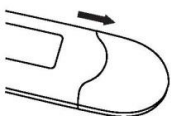
Figure 3

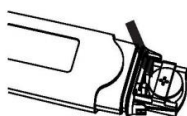
Figure 4

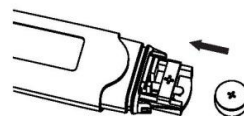
Figure 5

GUARANTEE AND CUSTOMER SERVICE

Before delivery our devices are subjected to rigorous quality control. If, despite all care, damage has occurred during production or transportation, please return the device to your dealer. In addition to statutory legal rights, the purchaser has an option to claim under the terms of the following guarantee:

For the purchased device we provide 2 years guarantee, commencing from the day of sale. If you have a defective product, you can directly go back to the point of purchase.

Defects which arise due to improper handling of the device and malfunctions due to interventions and repairs by third parties or the fitting of non-original parts are not covered by this guarantee. Always keep your receipt, without the receipt you can't claim any form of warranty. Damage caused by not following the instruction manual, will lead to a void of warranty, if this results in consequential damages then we will not be liable. Neither can we hold responsible for material damage or personal injury caused by improper use if the instruction manual is not properly executed. Damage to accessories does not mean free replacement of the whole appliance. In such case please contact our service department. Broken glass or breakage of plastic parts is always subject to a charge. Defects to consumables or parts subjected to wearing, as well as cleaning, maintenance or the replacement of said parts are not covered by the warranty and are to be paid.

ENVIRONMENT FRIENDLY DISPOSAL

Recycling – European Directive 2012/19/EU

This marking indicates that this product should not be disposed with other household wastes. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

BATTERY MUST BE RECYCLED OR DISPOSED OF PROPERLY. DO NOT OPEN. DO NOT DISPOSE OR IN FIRE OR SHORT CIRCUIT.

Manufacturer: Joytech Healthcare Co., LTD

Address: No. 365, Wuzhou Road, Yuhang Economic Development Zone, Hangzhou, China 311100

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Producent: Joytech Healthcare Co., LTD